

euroimmun 

From Revvity

# EUROLabOffice 4.0




# Digital laboratory organisation and intelligent interconnectivity for your diagnostic laboratory

- Flexibility for your daily routine using EUROLabOffice 4.0 as central interface in your diagnostic laboratory
- Simple, intuitive operation with only a few clicks – user-friendly interface for each individual work station
- Faster result generation – optimised workflows with minimal hands-on time
- Convenient automation process – continuous exchange with the laboratory information system (LIS) and connected instruments
- Full focus on the patient – clear presentation of all results with a comprehensive patient history including all immunofluorescence images on one page

## Intelligent integration – secure, fast and traceable

- Continuous exchange between ELO 4.0 and LIS for data availability at all times
- Inter-connection of all work stations and automated laboratory systems, including optimised connection of random access instruments
- Central interface for all detection methods: IFA, ELISA, blot and ChLIA
- Optimised access to all data, anytime and from anywhere

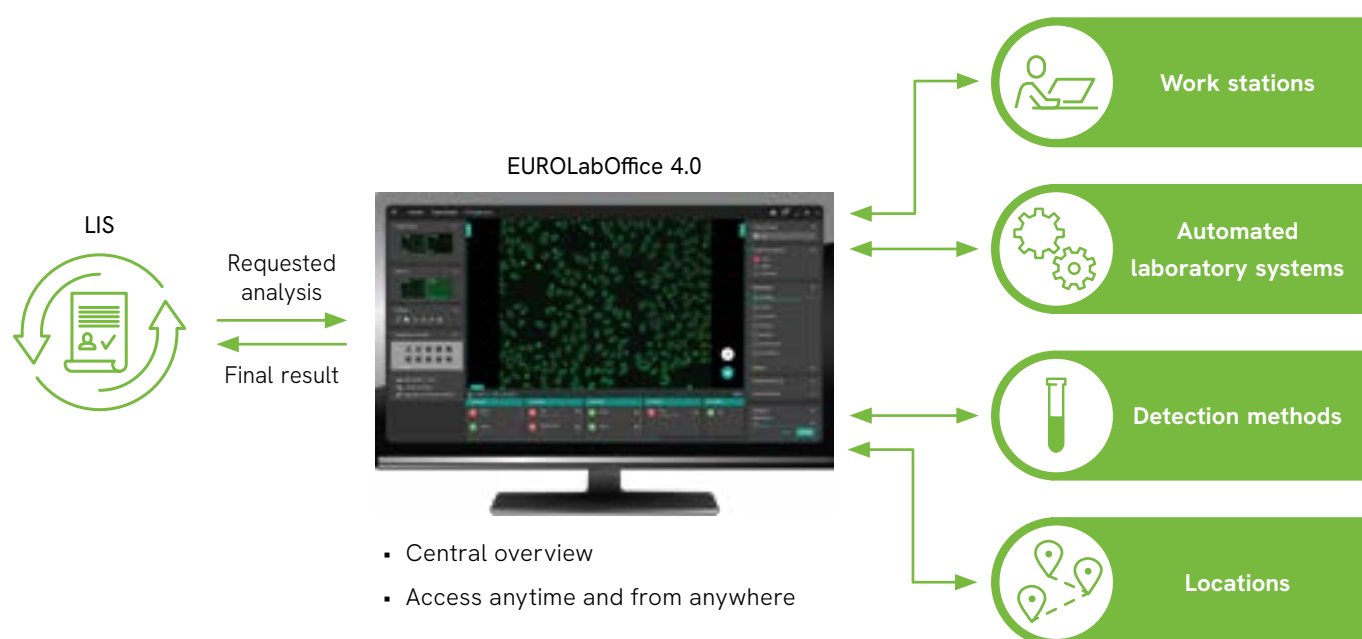

All in one software – valuable assistance from start to finish

1. Sample  
registration

2. Organisation



3. Test  
processing

4. Results

5. Archiving

## Flexible creation of optimised work protocols – for varying daily requirements in your laboratory routine

- Intuitive usage of the software thanks to the clearly structured user interface
- Automatic cost optimisation during protocol creation
- Automatic comparison of the samples registered with the laboratory instrument with open requests in the LIS
- Digital patient and result lists for simple laboratory organisation
- Possibility of spontaneous adaptation of sample, test or dilution choices via the live view and direct recalculation of the work protocol
- Flexible real-time access to all requests for optimal support in the laboratory routine

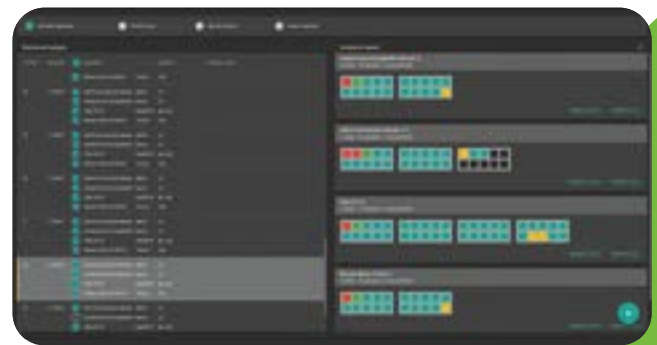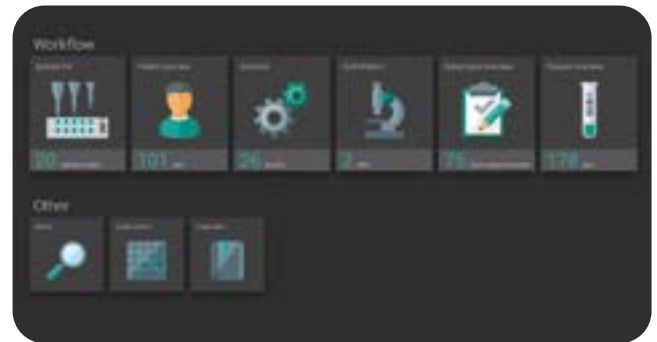

## User-friendly result display – fast and reliable result generation

- Helpful information on test processing and sample traceability available at all times
- Combined presentation of recent requests and all test results from all work stations on one page
- Access to the entire patient history, including IFA and blot images, directly during result entry
- Automatic result proposal for a growing number of substrates, encompassing positive/negative classification or pattern recognition as well as calculation of titers based on AI-enhanced software when using the UNIQO 160 or the automated microscopes EUROPattern Microscope Live and EUROPattern Microscope




- Result proposal



- Patient history with direct access to individual results and images

- Dilutions
- Substrates
- Process information
- Lot information

- Pattern and titer proposals
- Positive/negative classification

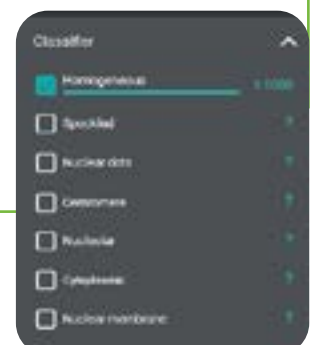



## The icing on the cake – maximum reliability of IFA evaluation with the image atlas



- Comprehensive IFA image atlas comprising ICAP\*-conformal reference images
- Automated selection of suitable reference images based on the software-suggested IFA pattern during use of the UNIQO 160 or the automated microscopes EUROPattern Microscope Live and EUROPattern Microscope
- Simple comparison of patterns on the evaluated image with the corresponding reference image in a side-by-side view
- Flexible expansion of the image atlas by user-own IFA images and freely creatable keywords according to individual needs

\* International Consensus on ANA Patterns

### Full focus on the patient

To enable fast and reliable result interpretation, ELO 4.0 combines all the results on one page, allowing a comprehensive view of a patient's data. The result evaluation is additionally supported by automatically generated result proposals and the image atlas.

This will give you more efficiency and reduce test reruns, resulting in cost- and time-savings for you and the patients.

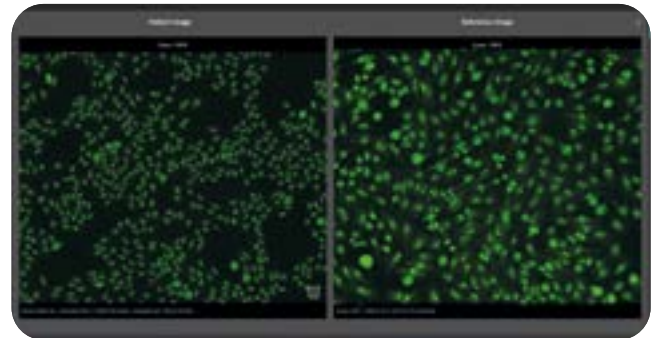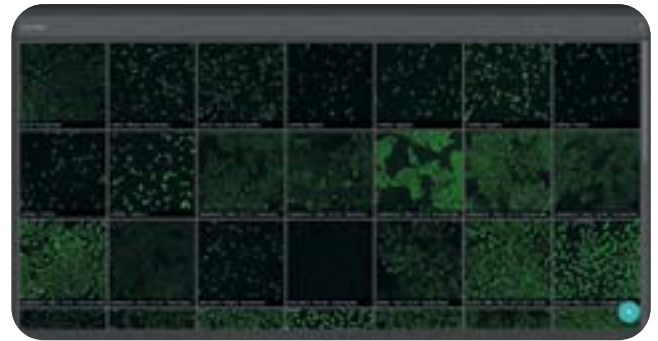

Comparison of patterns

## Instrument and software connections – a variety of solutions for every laboratory

- Bidirectional connection to fully automated IFA processing instruments, fully automated ELISA processors, automated microscopes as well as the fully automated random access systems (see Figure)
- Bidirectional connection to the automated blot processors EUROBlotOne and EUROBlotMaster (via EUROLineScan)
- Connection of devices from other manufacturers available upon request
- Standard interfaces with LIS connection (HL7, ASTM, GDT, XML etc.)



ELO  
4.0

EUROLabWorkstation IFA



EUROPattern Microscope



IF Sprinter



UNIQO 160



EUROPattern Microscope Live



Sprinter XL



EUROLabWorkstation ELISA



EUROArrayProcessor



EUROBlotOne



EUROIMMUN Analyzer I



IDS-i10\*\*



EUROBlotMaster 44



EUROIMMUN Analyzer I-2P



IDS-iSYS\*\*



EUROLineScan



\*\* Product manufactured by Immunodiagnostic Systems Ltd. (IDS). IDS-i10 and IDS-iSYS are trademarks of Immunodiagnostic Systems Ltd. (IDS), available to order from your local Euroimmun contact.

## Order information

Software	Order number
EUROLabOffice 4.0	YG 0250-0101-7

We are happy to configure EUROLabOffice 4.0  
according to your individual work needs.

Find out more at about EUROLabOffice 4.0  
[www.euroimmun.com](http://www.euroimmun.com)  
or contact us directly at  
[www.euroimmun.com/contact](http://www.euroimmun.com/contact)

